


### DYNA-FORM MERCURY STATIC MATTRESS



FEATURES	BENEFITS
Castellated foam insert (with tear duct cut)	Smaller castellated cells for heels and larger cells sacrally
Finest quality CME foam base and side walls	Prevents patient bottoming out and aids patient transfer
Multi-stretch, vapour permeable fabric (impregnated with an anti fungal agent)	Reduces shear & friction forces, whilst assisting infection control policies
Fully high frequency welded cover	Prevents fluid ingress & aids decontamination processes (can be laundered up to 95°C)
Toughened PU / non-slip base	Provides a robust protective surface against hydrolysis / strikethrough
Fully replaceable parts	Ensures great value for money

254KG SWL

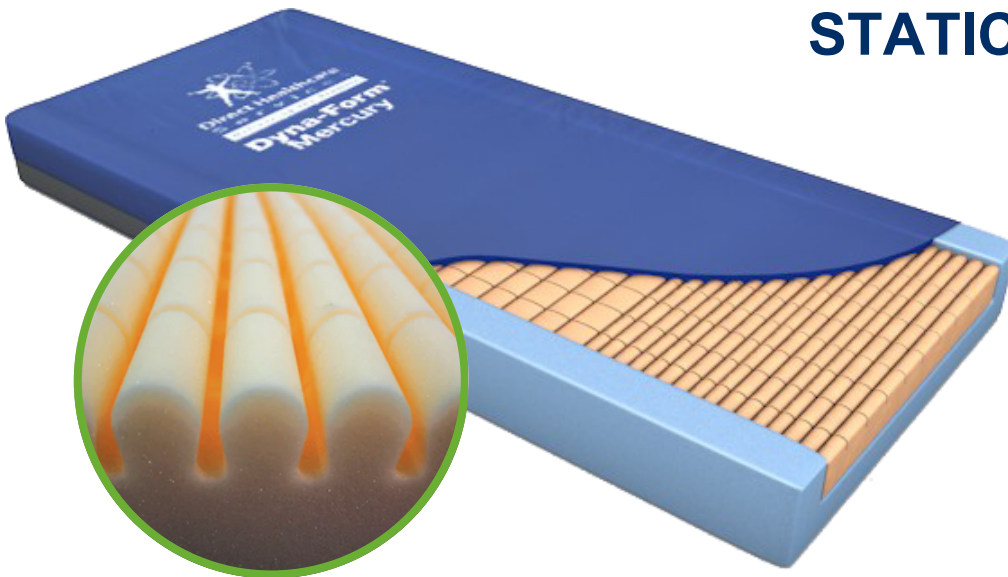

HIGH RISK



SUITABLE FOR STAGE 2 MANAGEMENT

2

8 YEAR FOAM WARRANTY



**CONTACT US**


The Mercury mattress is specifically designed to be used with patients considered to be at “High Risk” of pressure ulcer development. The extra rigidity of the Side Wall construction enables primary stability for patient transfer purposes. The “castellated cut” foam allows this product to provide excellent pressure re-distribution properties and ensure superior patient comfort and product longevity.

The Mercury is completed with a durable multi-stretch, vapour-permeable cover and a toughened polyurethane base. It provides an important contribution toward tissue viability standards of a modern hospital, nursing or home care environment.

Available in Single and King Single.

### KEY FEATURES

- Castellated foam insert: Tear duct cut
- Smaller cells for heels: Larger cells sacrally
- Quality CME foam base: Prevents patient bottoming out
- Side walls: Aids with patient transfer
- High frequency welded cover: Prevents fluid ingress
- Toughened PU base: Non-slip
- Multi-stretch, vapour permeable fabric: Reduces shear and friction forces

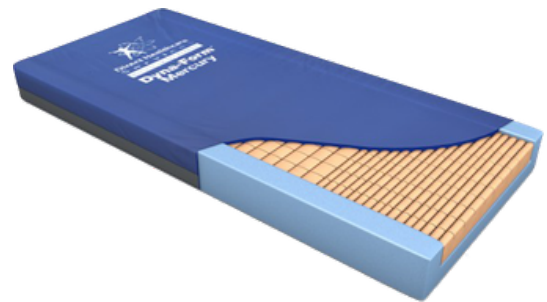

### SPECIFICATIONS

Model	Width	Length	Thickness	Product Weight	Warranty	SWL
Mercury Single	85.5cm	198cm	15cm	11kg	8 Year Foam / 4 Year Cover	254kg
King Single	107cm	198cm	15cm	13kg	8 Year Foam / 4 Year Cover	254kg

[CONTACT US](#)